



IconRehab

ICON *Portable Patient Hoists*



Optimised portable hoist design for safe patient transfers

An evolution in mobile hoists

The new **ICON Portable Hoist** family has been designed to offer a high feature-set of options that have been blended with new technologies to provide maximum comfort and security when lifting or transferring a patient to or from a bed, chair or even the floor.

Based on the user feedback, the new range has been customer insight led, encompassing everything that was valued in a mobile hoist whilst bringing in new features to further enhance the client and carer experience.

Impressive lifting range

The **ICON Portable Hoist** range encompasses an ergonomic lift that allows the client to be easily rotated 360°, making positioning much easier for the carer and reducing anxiety for the client.

Maximised space

The high boom ensures there is space for the client, even when the hoist is in its highest position, and the optimised space in front of the actuator reduces the risk of injury to the client's knees.



Easy to handle, transport and store

The new **ICON Portable Hoist** family can be easily folded and unfolded without the need for tools. They take up minimal space when stored and are easy to push and transport. If required, the mast can be very easily detached from the base (no tools required).

The **ICON Compact Portable Hoist** has a low weight of 31 kg making it even easier to transport and store.



NEW

**EXCLUSIVE
FEATURES**

**Curved Design
of Boom and
Spreader Bar**

- ▶ Distinct design allows carers to maintain eye contact with clients throughout a transfer.

**Ergonomic
Push Bar**

- ▶ Ergonomic design of the push bar makes handling and manoeuvring easier

**Battery Status
Indicator**

- ▶ Delivers peace of mind regarding the remaining lifts that are achievable

**Curved Leg
Design**

- ▶ Provides a cocoon like feeling and enables easy movement around objects.



Key features at a glance



◀ Hook Design

- Securely attach slings with ease, thanks to a new intuitive hook design
- Special ergonomic shape to attach sling from above or from below
- Flexible and outward pointing lip enables easy sling detachment



◀ Ergonomic Push Bar

- Ergonomic design of the push bar makes handling and manoeuvring easier
- Increased depth of handles allows for comfortable arm extension during transportation
- Designed with carers in mind



◀ Curved boom and spreader bar

- Distinct design of the boom and spreader bar allows carers to maintain eye contact with clients throughout a transfer
- Uninterrupted field of vision helps to alleviate any concerns for the client
- Enables 360° rotation of the spreader bar at the highest position to protect the leg from coming into contact with the mast



◀ Status indicator on control box

- Battery gauge indicator included as standard on all models
- Delivers peace of mind regarding the remaining lifts achievable
- Battery can be easily detached for efficient off-board charging
- *ICON PLUS* has Jumbo care control box as standard**

Introducing the **ICON** Portable Hoist Family

The mobile hoist range.
Designed with clients and carers in mind.



ICON
PORTABLE

▲ **The new standard in client transfer**

A redesigned version of market leading mobile hoists, designed to support safe transfers in both domestic and residential care.



ICON
PORTABLE COMPACT

▲ **Smaller footprint for easier manoeuvrability**

A smaller size allows the ICON PORTABLE COMPACT to easily operate in environments where space is limited.



ICON
PORTABLE PLUS

▲ **An advanced lifter with increased SWL**

Our plus sized model offers all the benefits of the ICON Portable with a higher SWL at 210kg, electrical legs and covered actuators as standard.

Lifting range	▶ 415 - 1660 mm	▶ 530 - 1590 mm	▶ 415 - 1660 mm
Safe working load	▶ 180 kg	▶ 150 kg	▶ 210 kg
Turning diameter	▶ 1400 mm	▶ 1070 mm	▶ 1430 mm


Available options and accessories:

- ▶ Electrical or mechanical leg spread (available on 180kg model only). Electrical is standard on **PLUS** model.
- ▶ Selection of 2 point or 4 point spreader bars in various widths.
- ▶ Extensive range of durable slings available in various sizes, styles and fabrics
- ▶ Lever for mechanical leg spread Homecare scale
- ▶ Leg protectors



Technical data ▼

For more comprehensive pre-sales information about this product, including the product's user manual, please visit www.iconrehab.co.uk.

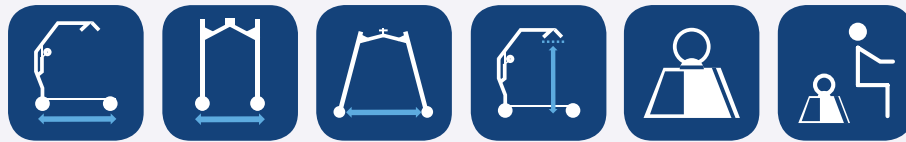

Base length ▼ Base width ▼ Inside width legs open ▼ Lifting range ▼ Product weight ▼ Safe working load ▼

Model	Base length	Base width	Inside width legs open	Lifting range	Product weight	Safe working load
ICON PORTABLE	1220 mm	640 mm	1020 mm	415 - 1660 mm	40 kg	180 kg
ICON PORTABLE COMPACT	1080 mm	520 mm	815 mm	530 - 1590 mm	31 kg	150 kg
ICON PORTABLE PLUS	1280 mm	640 mm	1020 mm	415 - 1660 mm	41 kg	210 kg



Height of legs ▼ Ground clearance ▼ Height when folded ▼ Turning diameter ▼ *Working ability ▼

Model	Height of legs	Ground clearance	Height when folded	Turning diameter	*Working ability
ICON PORTABLE	115 mm	35 mm	455 mm	1400 mm	40 cycles
ICON PORTABLE COMPACT	115 mm	35 mm	455 mm	1070 mm	40 cycles
ICON PORTABLE PLUS	115 mm	35 mm	470 mm	1430 mm	40 cycles


All measurements are taken with 100 mm castors and a 2 point spreader bar.
For 75 mm castors, reduce the height measurements by 15 mm.

* With battery charge between 100% and 50%
** Includes diagnostic functionality

Frame colour ▼



Light Grey:
mast, boom and spreader bar
RAL: 7035

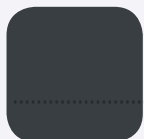

Anthracite Grey:
base and push bar
RAL: 7016

Regulatory data ▼



Please note that the model shown on the front cover is the 180kg model with electronically operated legs.

Warranty ▼

1 year warranty:

Frame, Actuator,
Control box, Hand
control, Battery

Please note: colours may vary slightly from those displayed above.